

# High Performance Cooled CCD Camera System ALTA U1109




The Alta U1109 has a back-illuminated full frame spectroscopy format CCD with exceptionally high quantum efficiency across the visible and in the UV. The 12 micron pixels allow high spectral resolution.

## Imaging Area of CCD



- Spectroscopy
- TDI Imaging


- 2048 x 506 array, 12 x 12 micron pixels
- 500 kHz 16-bit digitization
- 32Mbyte camera memory
- USB 2.0 interface: no plug in cards or external controllers
- Programmable, intelligent cooling to 55°C below ambient
- Binning up to 8 Horizontal x up to 506 Vertical
- Subarray readout and fast sequencing modes
- Precision time delayed integration (TDI) and kinetics mode readout
- Programmable fan speed for low / zero vibration
- Two serial port outputs for control of peripheral devices
- General purpose programmable I/O port
- External triggering and strobe controls
- ActiveX drivers and MaxIm DL/CCD software included with every system
- Field upgradeable firmware
- Fused silica windows
- Runs from single 12V supply with input voltage monitor
- Compact enclosure
- Programmable status indicators

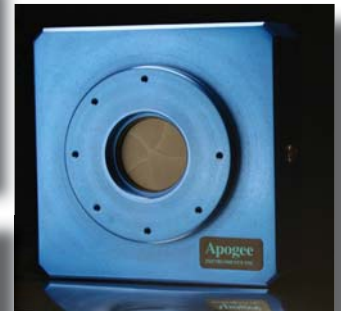

## CCD SPECIFICATIONS

CCD	Hamamatsu S10140-1109
Array Size (pixels)	2048 x 506
Pixel Size	12 x 12 microns
Imaging Area	24.6 x 6.1 mm (149.2 mm <sup>2</sup> )
Imaging Diagonal	25.3 mm
Video Imager Size	1.6"
Linear Full Well (typical)	75K electrons
Dynamic Range	81 dB
QE at 400 nm	58%
Peak QE (600-650 nm)	89%
Anti-blooming	none

For complete CCD specifications, including cosmetic grading, see data sheet from manufacturer.

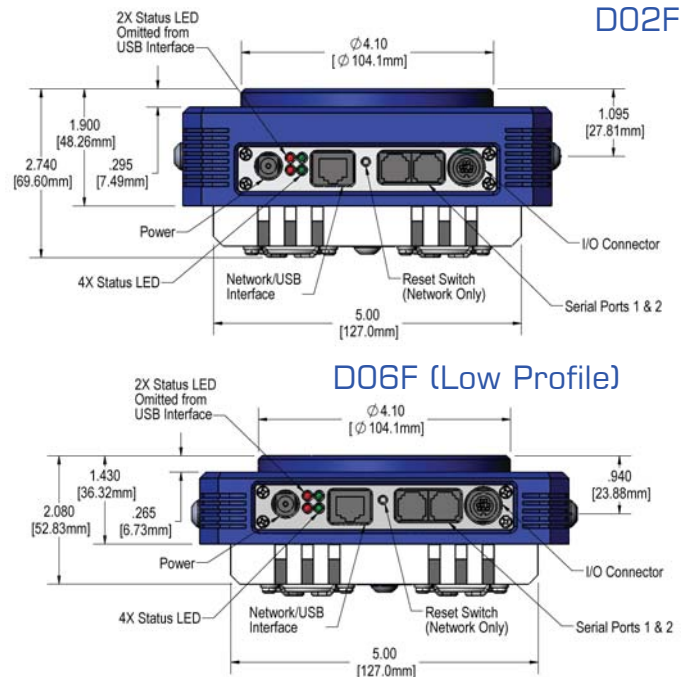

©2008 Apogee Instruments Inc. Alta is a registered trademark of Apogee Instruments Inc. Specifications subject to change without notice.




1020 Sundown Way, Ste 150  
Roseville CA 95661 USA  
tel 916-218-7450  
fax 916-218-7451  
<http://www.ccd.com>



PC Interface	USB 2.0
Max. Cable Length	5 meters between hubs; 5 hubs maximum (max. total of 30m)
Digital Resolution	16 bits at 500 kHz
System Noise (typical)	15 e <sup>-</sup> RMS at 500 kHz
Pixel Binning	1x1 to 8 x 506 on-chip
Exposure Time	30 milliseconds to 183 minutes (2.56 microsecond increments)
Image Sequencing	1 to 65535 image sequences under software control
Frame Sizes	Full frame, subframe, focus mode
Cooling (typical)	Thermoelectric cooler with forced air. Maximum cooling 50°C below ambient temperature
Dark Current (typical)	0.2 e <sup>-</sup> /pixel/sec (-30°C)
Temperature Stability	± 0.1°C
Camera Head Size	D2. Low profile: D6. Aluminum, hard blue anodized. 6" x 6" x 2.5" (15 x 15 x 6.35 cm) Weight: 3.1 lb. (1.4 kg)
Mounting	3.5" bolt circle. 2" 24 tpi thread. Optional Nikon F-mount or Canon FD mount.
Back Focal Distance	Standard: 1.025" (2.60cm). Low profile: 0.460" (1.17 cm). [optical]
Operating Environment	-22° to 27°C. Relative humidity: 10 to 90% non-condensing.
Cable Length	Standard: 15 ft (4.5m)
Power	40W maximum power with shutter open and cooling maximum. AC/DC "brick" supply with int'l AC input plug (100-240V, 50-60 Hz). Alternate 12V input from user's source.
Shutter	Standard: Melles Griot 43mm. Low profile: no shutter.
Remote Triggering	LVTTL input allows exposure to start within 25 microseconds of rising edge of trigger

### CCD SENSITIVITY

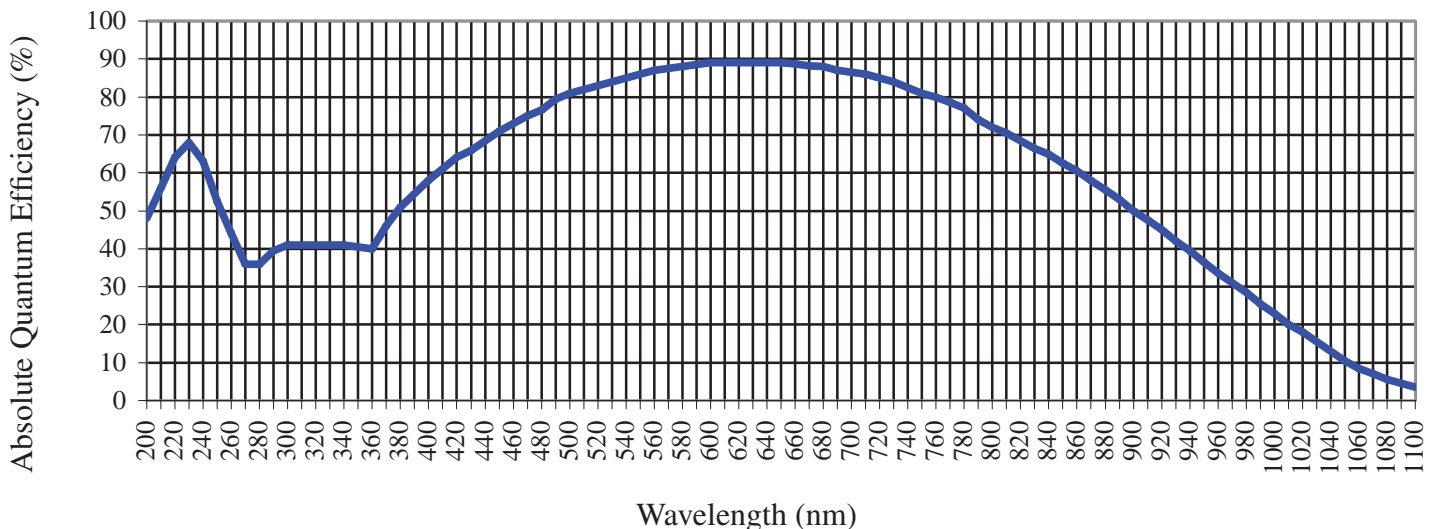

©2008 Apogee Instruments Inc. Alta is a registered trademark of Apogee Instruments Inc. Specifications subject to change without notice.